

Cassette holders**9890 010 8034 .****9890 010 8140 .****FILING INSTRUCTIONS**

File this documentation in binder:

SUBSYSTEM manual bucky DIAGNOST VE2 / VT2

UNIT manual bucky DIAGNOST VE / VT



PHILIPS

Philips Medical Systems Development and Manufacturing Centre

SERVICE MANUAL
704
UNIT

Cassette holders

9890 010 8034 .

9890 010 8140 .



DMC Hamburg

Printed in Hamburg, Federal Republic of Germany

© 2000 Philips Medical Systems
ALL RIGHTS RESERVED

Manual Order No. 4512 984 23621

Cassette holder

Cassette_holder_titel_000

(00.0)

SERVICE MANUAL - UNIT**Cassette holder**

Author: V. Neumann

Techn. No: 9890 010 8034 .
9890 010 8140 .

In case there are any questions concerning this manual,
please send this LOPAD via fax to 49/(0)40/5078 2481

File: Cassette holder (00.0)

List of pages and drawings (LOPAD)**Manual Order No: 4512 984 23621**

0.5	(00.0)	223 mm	(Rosa Karton)
1	(00.0)		
2	(00.0)		
3 . . .4	(00.0)		

1. Cassette holder


- Applies to bucky Diagnost VE / VT / VE2 / VT2
- Time to complete 1 engineer – 0.5 hours Tools : Standard set of tools
- Procedure:
- Determine the working position. Assure that both brackets "A" fits behind the profile rail. Correction: Readjust the brackets in horizontal direction.

Note:

The mark on extension rails are "setting marks" for bucky VE/VT front plate. Centre of cassette holder must be in centre of bucky wall stand.


- Support angle plate:

Support angle plate (B) grabs around the front plate of the wall stand. Correction: Readjust the support plate in vertical direction.

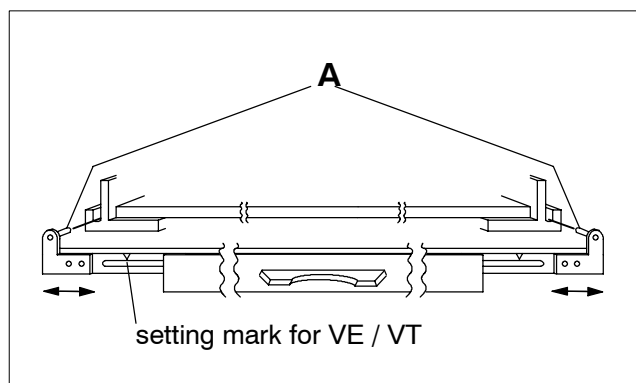


2. Function check:

- The cassette holder is mechanically adjusted if:
- Both brackets "A" are engaged in the profile rail
- Support angle plate "B" grabs around front plate
- The suckers "C" are engaged and the cassette holder is "tide".
- Remove the cassette holder: press both degassing studs "D" to loose the suckers "C". Tilt and lift the cassette holder.

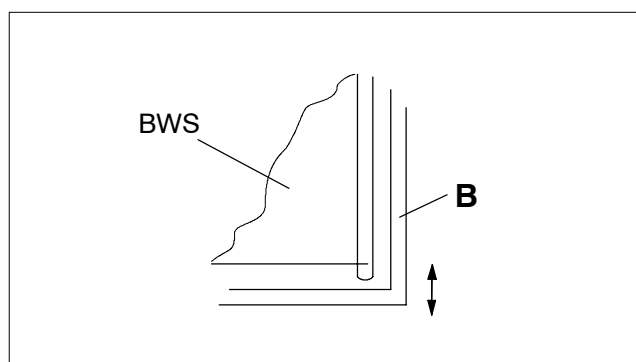


3. Exchange of a defect sucker / degassing stud:

Squeeze sucker from edge towards centre and turn around whilst pulling it off. Now the degassing stud can be taken out. The installation has be done in reverse order.

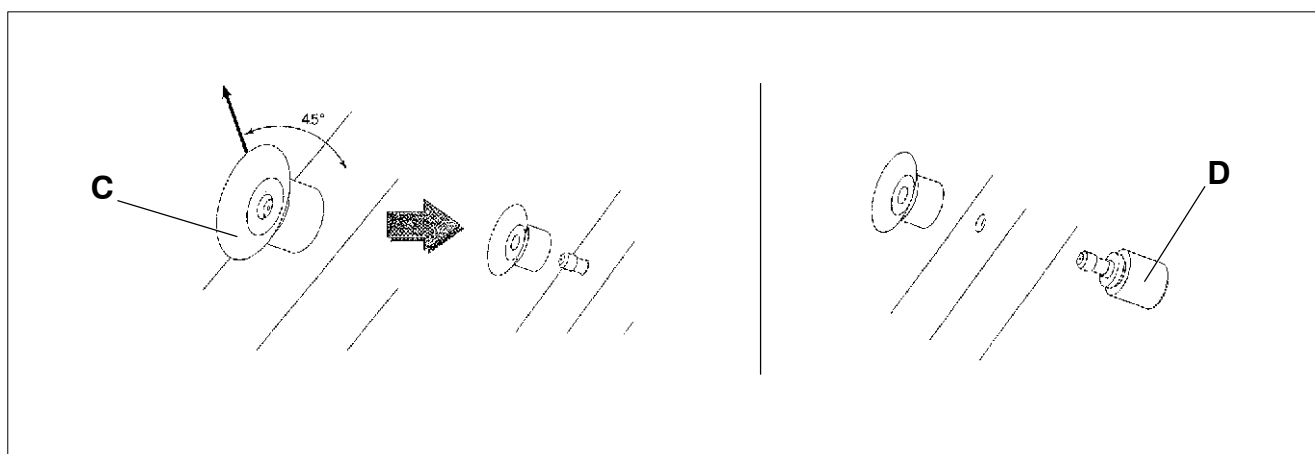


4. Labeling & Spare part list

 PHILIPS	Philips Medical Systems Development and Manufacturing Centre Röntgenstrasse 24 D-22335 Hamburg/Germany
type 9890 010 80342 s/n 00.00.123 cassetteholder	 0123

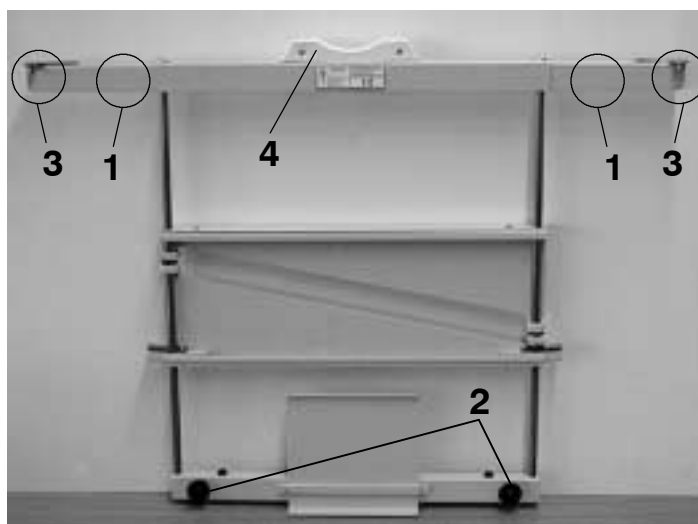

	9890 010 80342 Cassette holder normal
1	4512 520 39931 Support extension left / right
2	4512 520 39941 degassing stud assembly
3	4512 520 39961 Support angle left / right
4	4512 520 39921 Chin support assembly
5	4512 520 39951 Mounting kit complete (screws & nuts)
	4512 520 39991 Type lable
	4512 984 23621 Service manual

 PHILIPS	Philips Medical Systems Development and Manufacturing Centre Röntgenstrasse 24 D-22335 Hamburg/Germany
type 9890 010 81401 s/n 00.00.123 cassetteholder	 0123

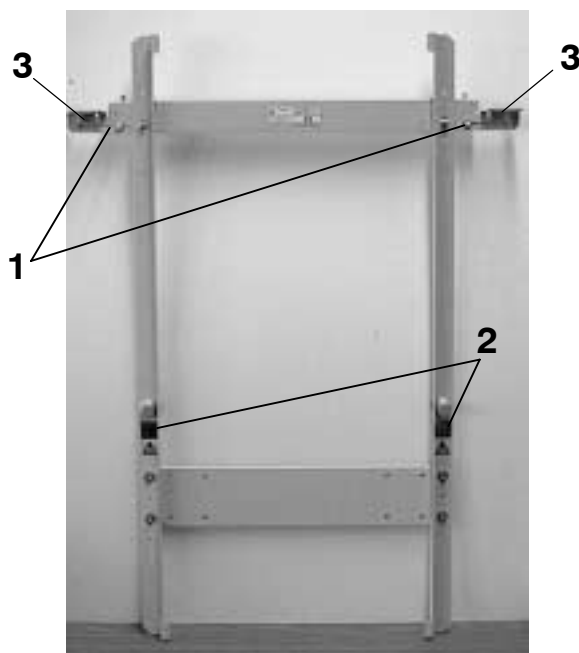

	9890 010 81401 Cassette holder long
1	4512 520 39911 Support extension left / right
2	4512 520 39941 degassing stud assembly
3	4512 520 39961 Support angle left / right
4	4512 520 39951 Mounting kit complete (screws & nuts)
	4512 520 39981 Type lable
	4512 984 23621 Service manual